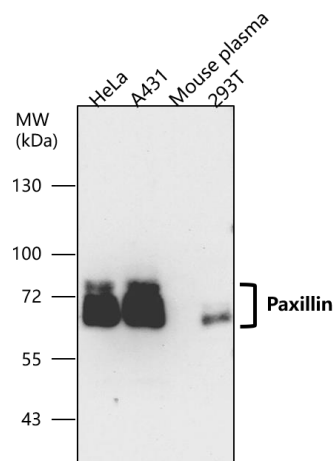


| | | |
|---|---|---|
| Datasheet for IR60-205 anti-Paxillin antibody | | Package size: 25, 100 µl Store at: -20°C MW (kDa): 68 |
| Overview | | |
| Product Name | Paxillin antibody | |
| Product Number | IR60-205 | |
| Gene Description | paxillin | |
| Clonality | Polyclonal | |
| Host | Rabbit | |
| Species Reactivity | Human , mouse, rat | |
| Recommended Applications Dilutions | Western Blot 1:1000 – 1:3000 Immunofluorescence 1:200 – 1:500 Immunohistochemistry (Paraffin) 1:300 – 1:400 | |
| Storage Buffer | 100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative | |
| Concentration | 1.01 mg/ml | |
| Purity | Affinity column purified | |
| Storage | Store at +4°C for short term storage. Long time storage is recommended at -20°C | |
| Notes | Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user. | |
| Image | | |
|  | | IR60-205 anti-Paxillin antibody WB image All lanes : Anti-Paxillin antibody at 1/2000 dilution Lysates/proteins at 30 µg per lane This blot was produced using a 7.5% SDS-PAGE. The gel was run at 140V for 50 minutes before being transferred onto a Nitrocellulose membrane at 18V for 60 minutes. The membrane was then blocked for an hour before being incubated with IR60-205 overnight at 4°C. |

IREal Biotechnology Co., Ltd.


7F., No. 9, Sec. 3, Zhonghua Rd., North Dist., Hsinchu City 300, Taiwan (R.O.C.)

Website : www.irealbio.com

TEL : 03-526-0005

E-mail : irealbio@irealbio.com

Customer feedback Image

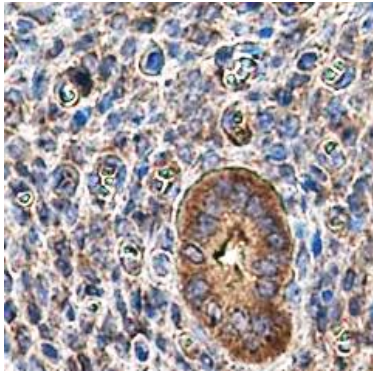

IR60-205 anti-Paxillin antibody WB image

Immunohistochemical of Paraffin sections.

Sample : mouse cervix.

Paxillin antibody (IR60-205) : 1-300.

IReal Biotechnology Co., Ltd.

7F., No. 9, Sec. 3, Zhonghua Rd., North Dist., Hsinchu City 300, Taiwan (R.O.C.)

Website : www.irealbio.com

TEL : 03-526-0005

E-mail : irealbio@irealbio.com

