

Catalog Number: IR131-528

Anti-XBP1 antibody

Package size: 25, 100 μ l

Store at: -20°C


MW (kDa): 60-80

Overview

Product Name	XBP1 antibody
Product Number	IR131-528
Gene Description	X-box binding protein 1
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, Mouse, Rat
Recommended Applications Dilutions	Western Blot 1:500-1500 Immunofluorescence 1:300 – 1:500
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.51 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

Image

IR131-528 anti-XBP1 antibody WB image

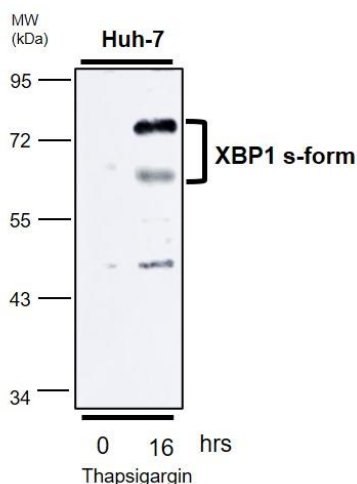

All lanes : Anti-XBP1 antibody - at 1/1000 dilution

Lysates/proteins at 50 μ g per lane

This blot was produced using a 10% SDS-PAGE. The gel was run at 140V for 50 minutes before being transferred onto a Nitrocellulose membrane at 18V for 60 minutes. The membrane was then blocked an hour before being incubated with IR131-528 overnight at 4°C.