


Catalog Number: IR250-898	Package size: 25, 100 µl
Anti-CD4 antibody	Store at: -20°C
	MW (kDa): 55-72


Overview

Product Name	CD4 antibody
Product Number	IR250-898
Gene Description	CD4 molecule
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:500 – 1:1000 Immunohistochemistry (Paraffin) 1:200 – 1:400 Flow Cytometry 1:100
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.51 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR250-898 anti-CD4 antibody WB image

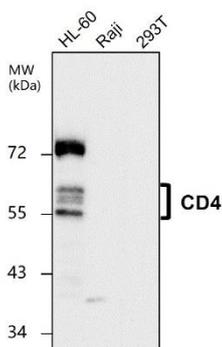

IR250-898 anti-CD4 antibody FACS image

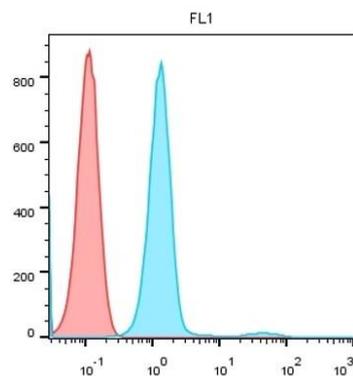

Anti-CD4 antibody at 1/700 dilution

Lysates at 50 µg per lane. This blot was produced using a 7.5% SDS-PAGE. NC membrane was then blocked for an hour before being incubated with IR250-898 overnight at 4°C.

Product: CD4 antibody

Catalog: IR250-898

Cell: Human peripheral blood lymphocytes

CD4 antibody: 1-100

Data from: POLOWANG BIOTECH INC

IReal Biotechnology Co., Ltd.