


Catalog Number: IR251-890
Anti-CD8 antibody

Package size: 25, 100 µl
 Store at: -20°C
 MW (kDa): 36


Overview

Product Name	CD8 antibody
Product Number	IR251-890
Gene Description	CD8a molecule
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:500 Immunohistochemistry (Paraffin) 1:200 – 1:400 Flow Cytometry 1:100
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.01 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR251-890 anti-CD8 antibody WB image

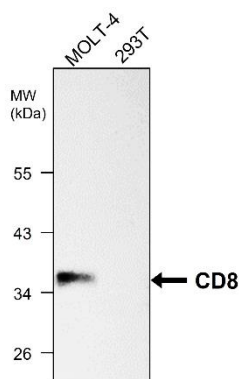

IR251-890 anti-CD8 antibody FACS image

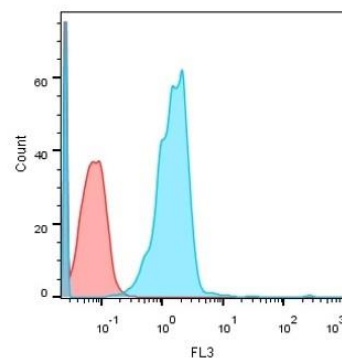

Anti-CD4 antibody at 1/500 dilution
 Lysates at 50 µg per lane. This blot was produced using a 10% SDS-PAGE. NC membrane was then blocked for an hour before being incubated with IR251-890 overnight at 4°C.

Product: CD8 antibody
 Catalog: IR251-890
 Cell: Human peripheral blood lymphocytes
 CD4 antibody: 1-100
Data from: POLOWANG BIOTECH INC

IReal Biotechnology Co., Ltd.