

Catalog Number: IR189-709


Anti-p38 MAPK antibody

Package size: 25, 100 μ l
 Store at: -20°C
 MW (kDa): 38-40


Overview

Product Name	p38 MAPK antibody
Product Number	IR189-709
Gene Description	mitogen-activated protein kinase 14
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat
Recommended Applications Dilutions	Western Blot 1:1000– 1:3000 Immunofluorescence 1:300 – 1:500 Immunohistochemistry (Paraffin) 1:100 – 1:200
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.5 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR189-709 anti-p38 MAPK antibody WB image

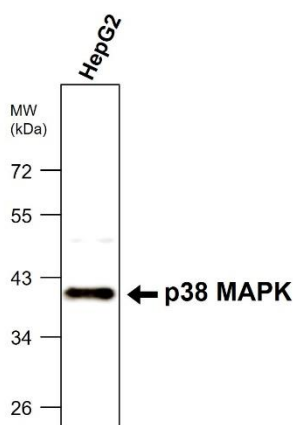

IR189-709 anti-p38 MAPK antibody IHC image

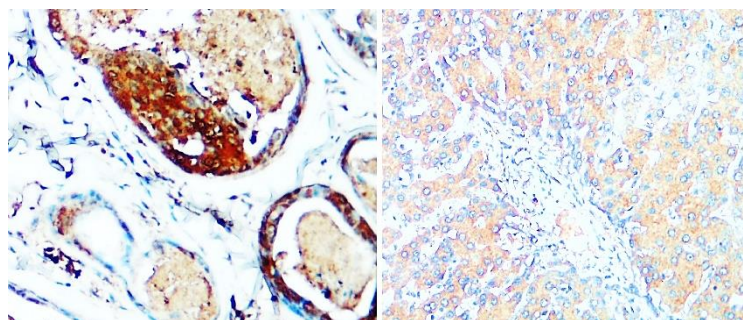

Anti-p38 MAPK antibody at 1/2000 dilution. Lysates at 30 μ g per lane. Nitrocellulose membrane was then blocked for an hour before being incubated with IR189-709 overnight at 4°C.

Immunohistochemical analysis of paraffin embedded Human cancer tissue labeling p38 MAPK with IR189-709 at 1/100.

IReal Biotechnology Co., Ltd.