

Catalog Number: IRT001

## Anti-HA tag antibody


Package size: 200, 100  $\mu$ l  
Store at: -20°C  
MW (kDa):

### Overview

Product Name	HA tag antibody
Product Number	IRT001
Clonality	Polyclonal
Host	Rabbit
Immunogen	Synthetic peptide (YPYDVPDYA)
Recommended Applications Dilutions	Western Blot 1:3000-1:7000 Immunofluorescence 1:500 – 1:1000
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.51 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

### Image

#### IRT001 anti HA-tag antibody WB image

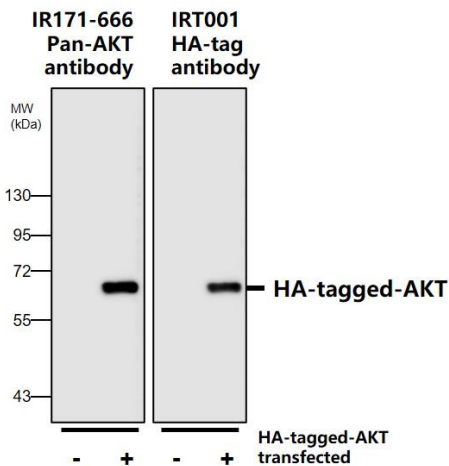

All lanes : Anti HA-tag antibody at 1/5000 dilution

Lysates/proteins at 30  $\mu$ g per lane

This blot was produced using a 7.5% SDS-PAGE. The gel was run at 140V for 50 minutes before being transferred onto a Nitrocellulose membrane at 18V for 60 minutes. The membrane was then blocked for an hour before being incubated with IRT001 overnight at 4°C.

IReal Biotechnology Co., Ltd.