

Catalog Number: IR291-959

Anti IRF-4 / MUM1 antibody

Package size: 25, 100 µl


Store at: -20°C

MW (kDa): 73


Overview

Product Name	IRF-4 / MUM1 antibody
Product Number	IR291-959
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, Mouse, Rat
Recommended Applications Dilutions	Western Blot 1:1000 - 1:3000 Immunofluorescence 1:300 - 1:600 Immunohistochemistry (Paraffin) 1:100 - 1:150
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.5 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR291-959 anti IRF-4 / MUM1 antibody WB image

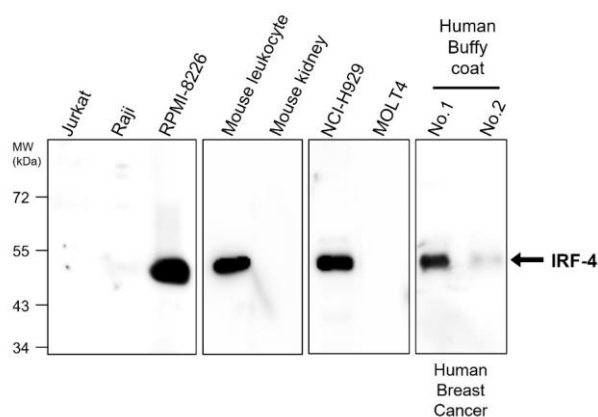

IR291-959 anti IRF-4 / MUM1 antibody IHC image

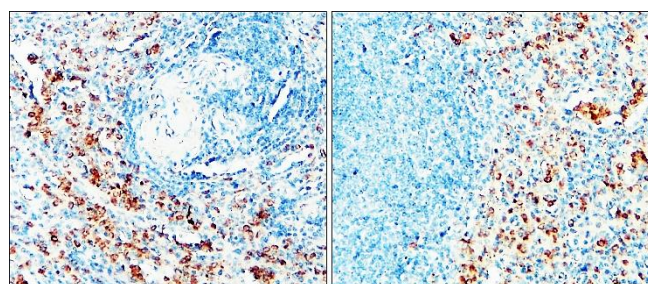

Anti IRF-4 / MUM1 antibody - at 1/2000 dilution

Lysates at 30 µg per lane. This blot was produced using a 7.5% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR291-959 overnight at 4°C.

Immunohistochemical analysis of paraffin embedded Human tonsil cancer tissue labeling IRF-4 / MUM1 with IR291-959 at 1/150.