

Catalog Number: IR023


Anti-Cyclophilin A / PPIA monoclonal antibody [23J]

Package size: 25, 100 µl
Store at: -20°C
MW (kDa): 17

Overview

Product Name	Cyclophilin A / PPIA monoclonal antibody [23J]
Product Number	IRM023
Gene Description	peptidylprolyl isomerase A (cyclophilin A)
Clonality	Monoclonal
Host	Mouse IgG1
Species Reactivity	Human, Mouse, Rat
Recommended Applications Dilutions	Western Blot 1:1000 Immunofluorescence 1:200 – 1:400
Storage Buffer	100mM Tris Glycine, 35% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1 mg/ml
Purity	Affinity column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IRM023 Cyclophilin A / PPIA monoclonal antibody [23J] WB image

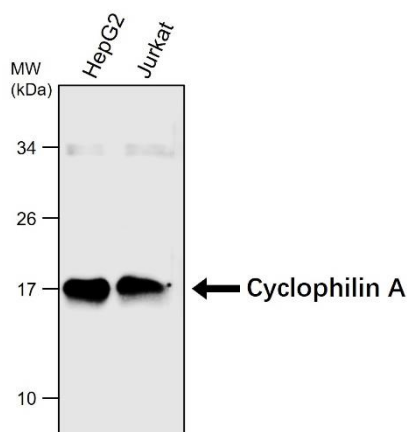

Anti-Cyclophilin A / PPIA, monoclonal antibody [23J] - at 1/1000 dilution

Lysates at 35 µg per lane.

This blot was produced using a 15% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IRM023 overnight at 4°C.

IReal Biotechnology Co., Ltd.