

Catalog Number: IRM214


Anti-Lactoferrin antibody [18-1]

Package size: 25, 100 μ l
 Store at: -20°C
 MW (kDa): 80-90

Overview

Product Name	Lactoferrin antibody [18-1]
Product Number	IRM214
Clonality	Monoclonal
Host	Mouse (mouse IgG1)
Species Reactivity	Human
Recommended Applications Dilutions	Western Blot 1:2000 – 1:5000 Immunofluorescence 1:500 – 1:1000
Storage Buffer	Phosphate-buffered solutions, pH 7.2, 50% Glycerol, 0.025% ProClin 300 was added as a preservative
Concentration	1mg/ml
Purity	Protein A column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

Image

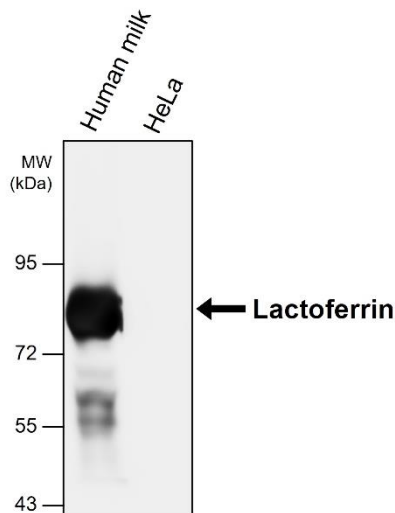

IRM214 Lactoferrin antibody [18-1] WB image

Anti-Lactoferrin antibody [18-1] at 1/2500 dilution

Lysates at 30 μ g per lane.
 This blot was produced using a 7.5% SDS-PAGE. The gel was run at 140V for 50 minutes before being transferred onto a Nitrocellulose membrane at 18V for 60 minutes. The membrane was then blocked for an hour before being incubated with IRM214 overnight at 4°C.