

Catalog Number: IRBS4032

Anti-Rat IgG (H&L) (GOAT) Antibody


Peroxidase Conjugated

Package size: 100, 200 µl
Store at: -20°C

Overview

Product Name	Anti-Rat IgG (H+L) (GOAT) Antibody Peroxidase Conjugated
Product Number	IRBS4032
Label	Peroxidase (Horseradish)
Clonality	Polyclonal
Host	Goat
Species Reactivity	Rat IgG
Recommended Applications Dilutions	Western Blot 1:5,000 - 1:10,000 ELISA 1:20,000 – 1:50,000 Other Assays User Optimized
Storage Buffer	0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2, 10 mg/mL Bovine Serum Albumin (BSA), 50% Glycerol, 0.025% ProClin 300 was added as a preservative, Do NOT add Sodium Azide!
Concentration	1 mg/ml by UV absorbance at 280 nm
Immunogen	Rat IgG whole molecule
Purity	Affinity column purified
Storage	Long time storage is recommended at -20°C
Application Note	<p>This anti-rat HRP secondary antibody is suitable for ELISA, Sandwich ELISA, titration assays, western blot, as well as other HRP antibody based assays.</p> <p>Specific conditions for reactivity and signal detection should be optimized by the end user. Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user. This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information.</p>

Image

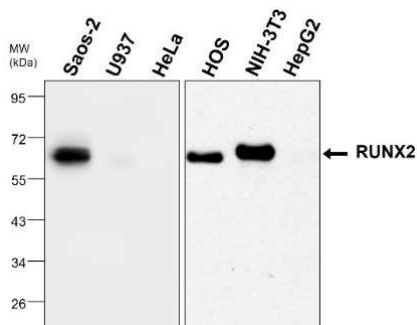

Western Blot of Peroxidase Conjugated Goat Anti-Rat IgG Antibody (**Catalog. IRBS4032**)

Lysates/proteins at 30 ug per lane
Primary antibody: IRR003 RUNX2, monoclonal antibody.

Secondary antibody: Peroxidase Conjugated Goat Anti- Rat IgG Antibody (**Catalog. IRBS4032**) at 1:6,000 for 60 min at RT.

IReal Biotechnology Co., Ltd.