

Catalog Number: IR43-128


Anti-alpha Synuclein / SNCA antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): 15


Overview

Product Name	alpha Synuclein / SNCA antibody
Product Number	IR43-128
Gene Description	synuclein, alpha
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat
Recommended Applications Dilutions	Western Blot 1:500 - 1:1000 Immunofluorescence 1:200 – 1:300 Immunohistochemistry (Frozen) 1:100 – 1:200
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR43-128 anti-alpha Synuclein antibody WB image

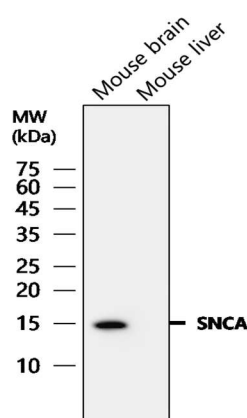

Customer feedback Image
IR43-128 anti-alpha Synuclein antibody IF image

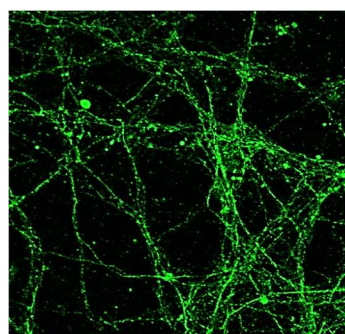

Anti-alpha Synuclein / SNCA antibody - 1/700 dilution

Lysates at 60 μ g per lane


This blot was produced using a 12% SDS-PAGE. The Nitrocellulose membrane was then blocked for an hour before being incubated with IR43-128 overnight at 4°C.

Immunofluorescent analysis.

Sample: primary cortical neurons
alpha Synuclein (IR43-128) : 1-400
Anti-rabbit 488: 1-500


Fixed: 4% paraformaldehyde at RT for 20 min.

Customer feedback Image
IR43-128 anti-alpha Synuclein antibody IHC image

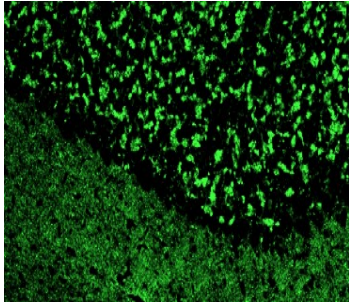

IR43-128 anti-alpha Synuclein antibody ICC/IF image

DBTRG-05MG

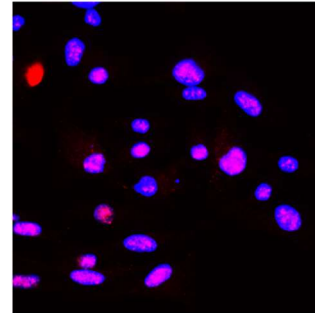

Immunohistochemical of frozen sections.

sample : mouse cerebellum.
alpha Synuclein (IR43-128) : 1-400
Anti-rabbit 488 : 1-500

Immunofluorescence: cells were fixed with 4% paraformaldehyde for 10 min at RT, permeabilized with 0.1% NP-40 for 10 min at RT then blocked with 5% BSA for 30 min at room temperature. Cells were stained with IR43-128 anti-alpha Synuclein antibody (red) at 1:200 and 4°C. DAPI (blue) was used as the nuclear counter stain.

iReal Biotechnology Co., Ltd.

Website : www.irealbio.com

TEL : 03-526-0005

E-mail : irealbio@irealbio.com

 **iReal Antibodies**



**ISO 13485:2016 Quality
Management System
Certified**

 www.irealbio.com

 03-5260005

 irealbio@irealbio.com