

Catalog Number: IRBS4012N

RealBlot anti-Native rabbit IgG antibody, HRP conjugated


Package size: 100 µl
Store at: -20°C

Overview

Product Name	RealBlot anti-Native rabbit IgG antibody, HRP conjugated
Product Number	IRBS4012N
Label	Peroxidase (Horseradish)
Clonality	Monoclonal
Host	Mouse
Species Reactivity	Native Rabbit IgG
Recommended Applications Dilutions	User Optimized Western Blot 1:3000 – 1:5000
Storage Buffer	100mM sodium phosphate(pH7.4), 40% Glycerol, 0.025% ProClin 300 was added as a preservative, Do NOT add Sodium Azide!
Concentration	1 mg/ml
Immunogen	Rabbit IgG whole molecule
Purity	Protein A column purified
Storage	Long time storage is recommended at -20°C

Image

RealBlot anti-Native rabbit IgG antibody, HRP conjugated

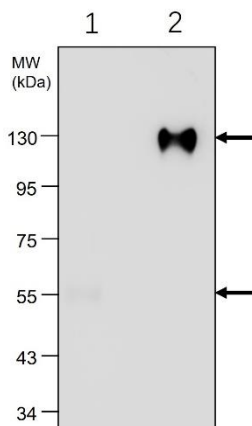

RealBlot anti-Native rabbit IgG antibody, HRP conjugated detects native-rabbit IgG by western blot.

Protein Loading :

1. 5ug Denature Rabbit IgG
2. 1ug Native Rabbit IgG

Protein was separated by 10% SDS-PAGE, and the membrane was blotted with RealBlot anti-Native rabbit IgG antibody, HRP conjugated (IRBS4012N) at a dilution of 1:3000.