


## Catalog Number: IR299-1 Anti Phospho-Troponin I (Cardiac) (Ser23/24) antibody

Package size: 25, 100  $\mu$ l  
Store at: -20°C  
MW (kDa): 28


### Overview

Product Name	Phospho-Troponin I (Cardiac) (Ser23/24) antibody
Product Number	IR299-1
Gene Description	troponin I, cardiac 3
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat
Recommended Applications Dilutions	Western Blot 1:1000 – 1:2000 Immunofluorescence 1:200 – 1:400
Storage Buffer	100mM Tris Glycine, 30% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	0.35 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

#### IR299-1 anti-Phospho-Troponin I (Cardiac) (Ser23/24) antibody WB image

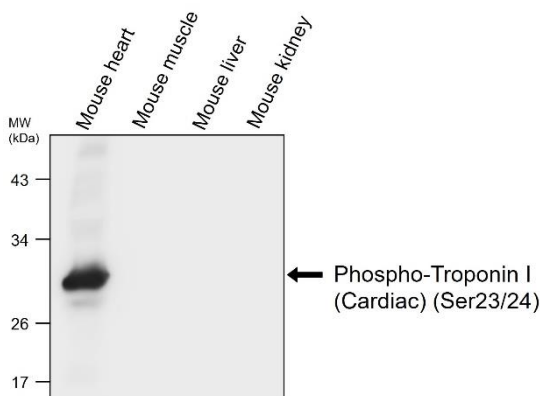

#### IR299-1 anti-Phospho-Troponin I (Cardiac) (Ser23/24) antibody WB image

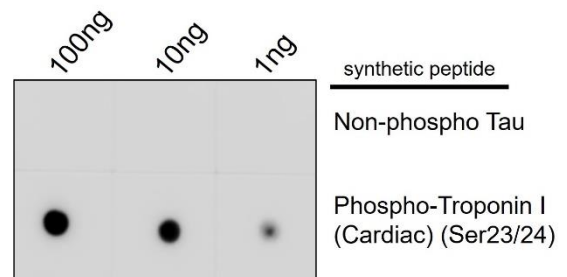

#### Anti-Phospho-Troponin I (Cardiac) (Ser23/24) antibody at 1/1000 dilution

Lysates/proteins at 20  $\mu$ g per lane  
This blot was produced using a 15% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR299-1 overnight at 4°C.

Dot blot for specificity and sensitivity of **Phospho-Troponin I (Cardiac) (Ser23/24) Antibody (IR299-1)**

IReal Biotechnology Co., Ltd.