

Product datasheet

WB

IF

Catalog Number: IRY25002

Anti-Prohibition / PHB IgY

Package size: 25, 100 µl


Store at: -20°C

MW (kDa): 29

Overview

Product Name	Prohibition / PHB IgY
Product Number	IRY25002
Clonality	Polyclonal IgY
Host	Chicken
Species Reactivity	Human, Mouse, Rat
Immunogen	Recombinant protein specific to human PHB1 protein
Recommended Applications Dilutions	Western Blot 1:5000 – 1:7000 Immunofluorescence 1:500 – 1:1000
Storage Buffer	Phosphate buffer, 30% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative.
Concentration	3 mg/ml
Purity	Eggs collected and isolated the IgY fraction (purity 90-95%).
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

Prohibition / PHB IgY WB image

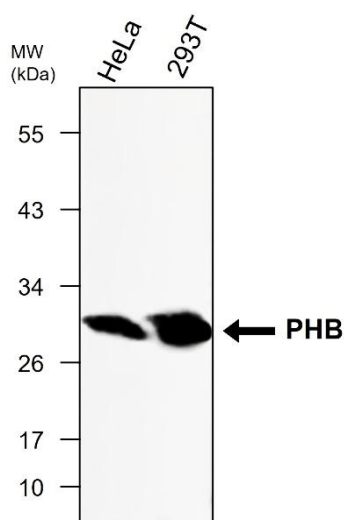

Anti-Prohibition / PHB IgY - at 1/5000 dilution
Lysates at 30 µg per lane.
This blot was produced using a 12% SDS-PAGE.
Nitrocellulose membrane was then blocked for an hour before being incubated with IRY25002 overnight at 4°C.

IReal Biotechnology Co., Ltd.

