

Product datasheet

WB

IF

Catalog Number: IRY25003

Anti-STIM1 IgY

Package size: 25, 100 µl


Store at: -20°C

MW (kDa): 85-90

Overview

Product Name	STIM1 IgY
Product Number	IRY25003
Clonality	Polyclonal IgY
Host	Chicken
Species Reactivity	Human, Mouse, Rat
Immunogen	Recombinant protein specific to human STIM1 protein
Recommended Applications Dilutions	Western Blot 1:5000 – 1:7000 Immunofluorescence 1:500 – 1:1000
Storage Buffer	Phosphate buffer, 30% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative.
Concentration	3 mg/ml
Purity	Eggs collected and isolated the IgY fraction (purity 90-95%).
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

STIM1 IgY WB image

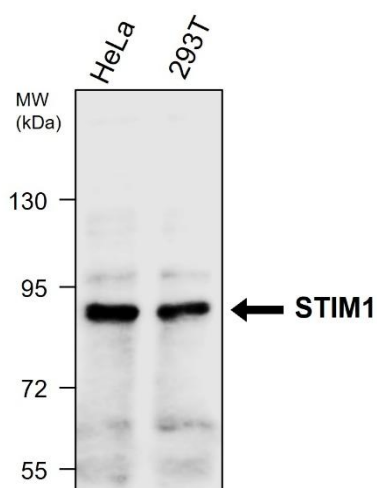

Anti-STIM1 IgY - at 1/5000 dilution

Lysates at 30 µg per lane.

This blot was produced using a 7.5% SDS-PAGE.

Nitrocellulose membrane was then blocked for an hour before being incubated with IRY25003 overnight at 4°C.

IReal Biotechnology Co., Ltd.

