

Catalog Number: IR60-205


Anti-Paxillin antibody

Package size: 25, 100 µl
Store at: -20°C
MW (kDa): 68


Overview

Product Name	Paxillin antibody
Product Number	IR60-205
Gene Description	paxillin
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat
Recommended Applications Dilutions	Western Blot 1:1000 – 1:3000 Immunofluorescence 1:200 – 1:500 Immunohistochemistry (Paraffin) 1:300 – 1:400
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	0.7 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR60-205 anti-Paxillin antibody WB image

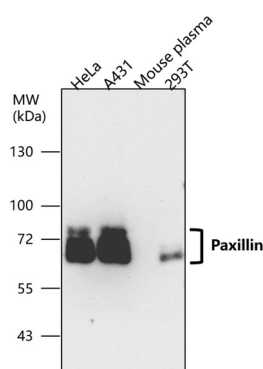

IR60-205 anti-Paxillin antibody IHC image

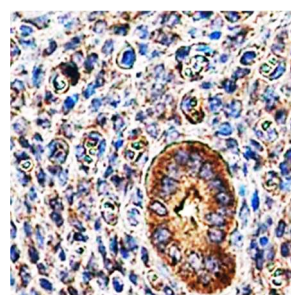

Anti-Paxillin antibody at 1/2000 dilution

Lysates/proteins at 30 µg per lane. This blot was produced using a 7.5% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR60-205 overnight at 4°C.

Immunohistochemical of Paraffin sections.
Sample: mouse cervix.
Paxillin antibody (IR60-205) : 1-300.

iReal Biotechnology Co., Ltd.