

Catalog Number: IR3-8


Anti-GAPDH antibody, HRP conjugated

Package size: 25, 100 μ l
 Store at: -20°C
 MW (kDa): 34

Overview

| | |
|-----------------------|---|
| Product Name | GAPDH antibody, HRP conjugated |
| Product Number | IR3-8-HRP |
| Conjugation | HRP conjugated |
| Clonality | Polyclonal |
| Host | Rabbit |
| Species Reactivity | Human, mouse, rat |
| Recommended Dilutions | Western Blot 1:5000 - 1:10000 |
| Storage Buffer | Phosphate buffered saline (pH7.4), 50% Glycerol. 0.025% ProClin300 was added as a preservative. |
| Concentration | 0.5 mg/ml |
| Form | Liquid |
| Purity | Affinity column purified |
| Storage | Long time storage is recommended at -20°C, Store in the dark. |
| Notes | Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user. |

IR3-8-HRP anti-GAPDH antibody, HRP conjugated WB image

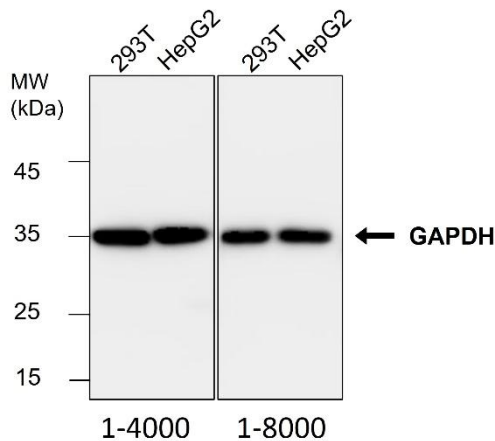

Lysates at 25 μ g per lane. This blot was produced using a 10% SDS-PAGE. Nitrocellulose membrane was blocked for an hour before being incubated with IR3-8-HRP one hour at room temperature. Antibody binding was visualized using ECL development solution. Exposure time: 10s

