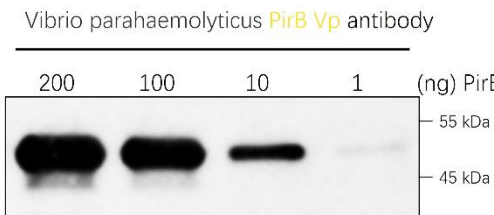
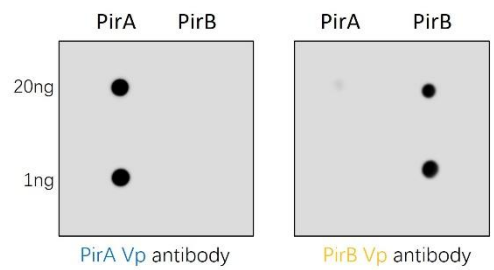


Catalog Number: IRB2511B	Package size: 25, 100µl Store at: -20°C
Anti-Vibrio parahaemolyticus PirB Vp antibody	

Overview

Product Name	Vibrio parahaemolyticus PirB Vp antibody
Product Number	IRB2511B
Immunogen	E. coli extract with Vibrio parahaemolyticus rPirB Vp expression.
Clonality	Polyclonal
Host	Rabbit
Recommended Dilutions	Western Blot 1:3000 – 1:6000 Optimal dilution to be determined by titration.
Storage Buffer	100mM Tris Glycine, 40% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1 mg/ml
Purity	Affinity column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

Image

IRB2511B anti-Vibrio parahaemolyticus PirB Vp antibody WB image	IRB2511B anti-Vibrio parahaemolyticus PirB Vp antibody WB image
<p style="text-align: center;">Vibrio parahaemolyticus PirB Vp antibody</p> 	
<p>Anti-Vibrio parahaemolyticus PirB Vp antibody at 1/5000 dilution</p> <p>This blot was produced using a 10% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IRB2511B overnight at 4°C.</p>	<p>Dot blot for specificity and sensitivity of Vibrio parahaemolyticus PirB Vp antibody (IRB2511B)</p>

IReal Biotechnology Co., Ltd.

