

WB

IF

IHC

Catalog Number: IR296-1

Anti-Collagen II / COL2A1 antibody

Package size: 25, 100 μ l


Store at: -20°C

MW (kDa): 180-200


Overview

Product Name	Collagen II / COL2A1 antibody
Product Number	IR296-1
Gene Description	collagen, type II
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat
Recommended Applications Dilutions	Western Blot 1:1000 Immunofluorescence 1:300 – 1:500 Immunohistochemistry (Paraffin) 1:100 – 1:200
Storage Buffer	100mM Tris Glycine, 30% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.36g/ml
Purity	Affinity column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR296-1 anti-Collagen II / COL2A1 antibody WB image

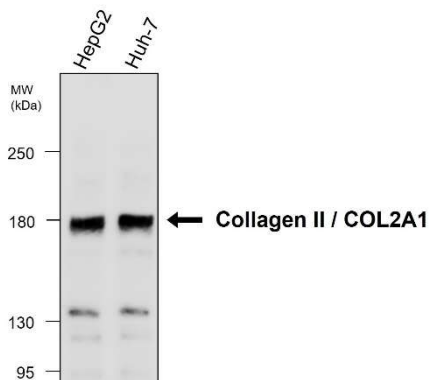

IR296-1 anti-Collagen II / COL2A1 antibody IHC image

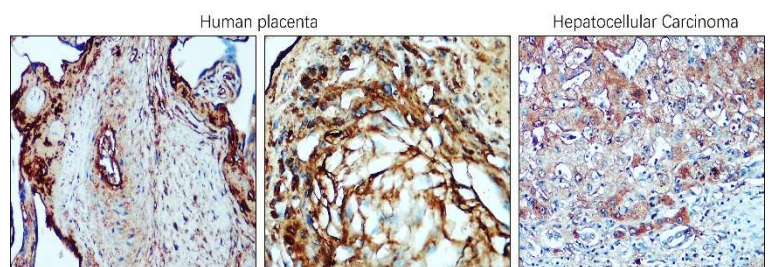

Collagen II / COL2A1 antibody - at 1/1000 dilution.

Lysates at 20 μ g per lane. This blot was produced using a 5% SDS-PAGE. Nitrocellulose membrane was then blocked before being incubated with IR296-1 overnight at 4°C.

Immunohistochemical analysis of paraffin embedded human tissue labeling **Collagen II / COL2A1** with IR296-1 at 1/100.

IReal Biotechnology Co., Ltd.

