

Catalog Number: IR359-1


Anti-Mucin 2 / MUC2 antibody

Package size: 25, 100 μ l
Store at: -20°C
MW (kDa): >350


Overview

Product Name	Mucin 2 / MUC2 antibody
Product Number	IR359-1
Gene Description	mucin 2, oligomeric mucus/gel-forming
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat
Recommended Applications Dilutions	Western Blot 1:1000 – 1:2000 Immunofluorescence 1:400 – 1:600 Immunohistochemistry 1:100 – 1:300
Storage Buffer	100mM Tris Glycine, 30% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.16 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR359-1 anti-Mucin 2 / MUC2 antibody WB image

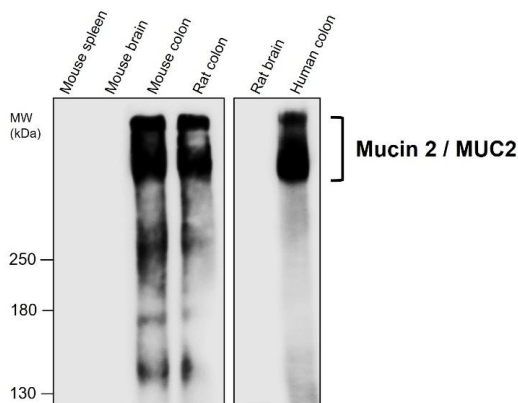

IR359-1 anti-Mucin 2 / MUC2 antibody IHC image

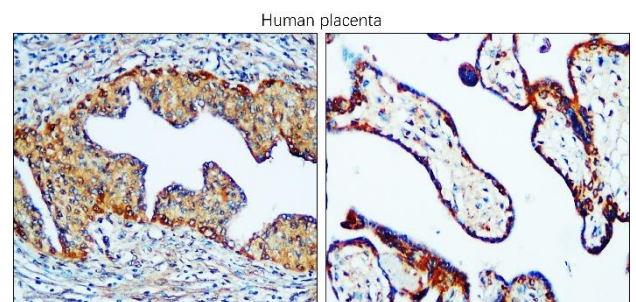

Anti-Mucin 2 / MUC2 antibody at 1/1000 dilution

Lysates/proteins at 40 μ g per lane
This blot was produced using a 5% SDS-PAGE. Nitrocellulose membrane was then blocked for an hour before being incubated with IR359-1 overnight at 4°C.

Immunohistochemical analysis of paraffin embedded human placenta labeling Mucin 2 / MUC2 with IR359-1 at 1/100.

IReal Biotechnology Co., Ltd.