

Catalog Number: IR282-B1

Anti CDH4 / R-cadherin antibody

Package size: 25, 100 µl


Store at: -20°C

MW (kDa): 110-130


Overview

Product Name	CDH4 / R-cadherin antibody
Product Number	IR282-B1
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat, bovine
Recommended Applications Dilutions	Western Blot 1:1000– 1:2000 Immunofluorescence 1:300 – 1:500 Immunohistochemistry (Paraffin) 1:100
Storage Buffer	100mM Tris Glycine, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.51 mg/ml
Purity	Affinity column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR282-B1 anti CDH4 / R-cadherin antibody WB image

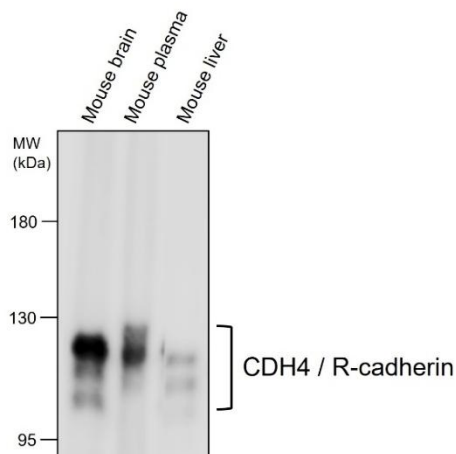

IR282-B1 anti CDH4 / R-cadherin antibody IHC image

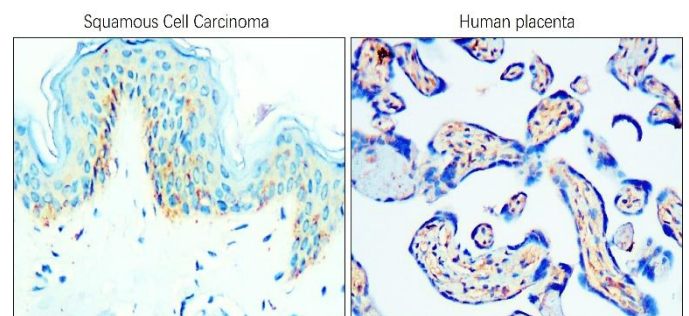

Anti CDH4 / R-cadherin antibody at 1/1000 dilution

Lysates at 30 µg per lane. This blot was produced using a 5% SDS-PAGE. The gel was run at 140V for 50 minutes before being transferred onto a Nitrocellulose membrane was then blocked for an hour before being incubated with IR282-B1 overnight at 4°C.

Immunohistochemical analysis of paraffin embedded human tissue labeling **CDH4 / R-cadherin** with IR282-B1 at 1/100.