

Catalog Number: IRM061

Anti-Iba1 monoclonal antibody [26F1]

Package size: 100 µl


Store at: -20°C

MW (kDa): 15~17


Overview

Product Name	Iba1 monoclonal antibody [26F1]
Product Number	IRM061
Clonality	Monoclonal
Host	Mouse (mouse IgG1)
Species Reactivity	Human, Mouse, Rat
Recommended Applications Dilutions	Western Blot 1:500 – 1:1000 Immunofluorescence 1:100 – 1:200 Immunohistochemistry (Paraffin) 1:100 – 1:200
Storage Buffer	Phosphate-buffered solution, pH 7.2, 40% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	1.2 mg/ml
Purity	Affinity column purified
Storage	Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IRM061 anti-Iba1 monoclonal antibody [26F1] IHC image

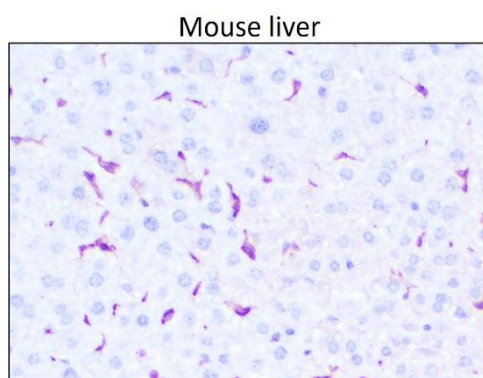

IRM061 anti-Iba1 monoclonal antibody [26F1] WB image

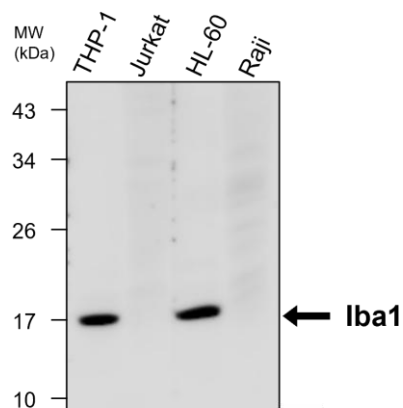

Iba1 monoclonal antibody [26F1] IRM061. Immunoperoxidase staining of formalin fixed, paraffin-embedded mouse liver tissue.

Iba1 antibody at 1/500 dilution

Lysates/ at 40 µg per lane. Nitrocellulose membrane was then blocked for an hour before being incubated with IRM061 overnight at 4°C.

IReal Biotechnology Co., Ltd.